

FLOTEC INC

7625 West New York Street
Indianapolis, Indiana 46214-4911

Phone: 317-273-6960
Fax 317-273-6979

Web Site: <http://www.floteco2.com>
E-mail: flotec@floteco2.com


Order desk: 800-401-1723
Order desk Fax 800-515-9254

Esophageal Detection Device (EDD)

Features

Easy to use

Positube attaches to the endotracheal tube post intubation. Resistance when drawing back on the syringe-type plunger indicates esophageal placement. When properly placed, the rigidity of the trachea allows air flow and creates no resistance when drawing back.

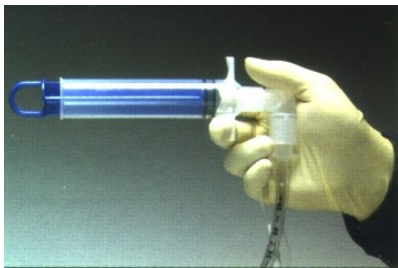

Ergonomic design

Pistol grip design with unique lateral syringe pull promotes better stability and enhanced tube positional control.

Positive Control

Patented design features finger grips on the proximal end ensuring stability while plunger is pulled. Unique 90 degree detachable elbow allows for lateral pull versus vertical preventing in-line tube movement or displacement. Fits any standard bag valve mask.

Large Pull Ring Accommodates gloved fingers easily

Modular Design

90 degree elbow attachment, while promoting better control and stability also prevents detachment of tube end adapter and allows for more convenient resuscitation bag attachment and positioning



Ordering Information	Description
PRA-PT	Esophageal Detection Device