

# FLOTEC

7625 West New York Street Indianapolis IN 46214-4911 [www.floteco2.com](http://www.floteco2.com)  
Phone: 317-273-6960 Fax: 317-273-6979 Order Desk: 800-401-1723  
e-mail: [flotec@floteco2.com](mailto:flotec@floteco2.com) Order Desk Fax: 1-800-515-9254

# Winner!

## EMS EXPO **TOP!** PRODUCT AWARD

# MASS CASUALTY MULTIPLE OXYGEN FLOW METER DISTRIBUTION ASSEMBLY

Now with Optional Swivel-T Connections



As shown above  
**Cat. # MCA-226P7-20FX0**

Flow Meter Chain Assembly with multiple (0-15) or (0-25) LPM adjustable flowmeters attached to flexible hose at 6 foot intervals using DISS-1240 Female Hand-Tight Connections.

Includes CGA-540 Pressure Reducing Regulator and 20 foot Oxygen Leader Supply Hose.

All-Metal-Constructed Flow Meters are extremely durable and are non-gravity sensitive.

Check-Valve Outlets can also be used to power CPAP Units and Resuscitators.


Overall extended length is 44ft. with 5 flowmeters, 50ft. with 6 flowmeters.

Flowmeter Drop Hoses available for Mobile Shelter / Hospital applications.


Customizable to your specific departmental requirements.

Connectable to any 50 psi oxygen gas source.


# MCA Series Multiple Oxygen Flow Meter Distribution Assemblies

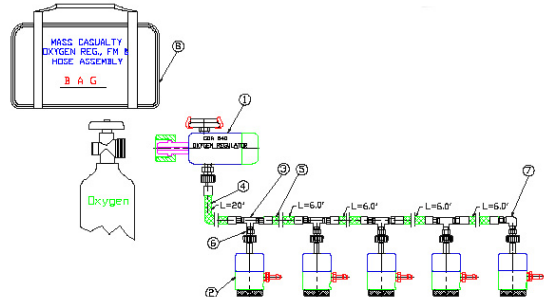


**MCA- 125P7-204X0**  
**MCA-125P7-20DX0 (with Swivel T's)**  
 5-Flow Meters (0-25 LPM) @ 6ft Intervals  
 CGA-540 Pressure Reducing Regulator  
 20ft Supply Feed Hose and Carry / Storage Bag

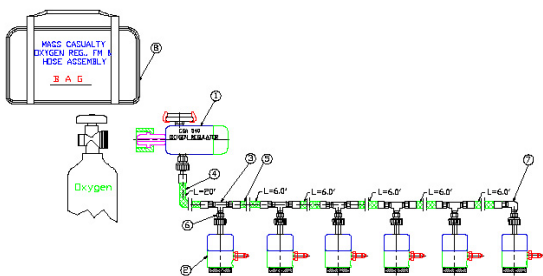

**MCA- 126P7-205X0**  
**MCA-126P7-20FX0 (with Swivel T's)**  
 6-Flow Meters (0-25 LPM) @ 6ft Intervals  
 CGA-540 Pressure Reducing Regulator  
 20 ft Supply Feed Hose and Carry / Storage Bag

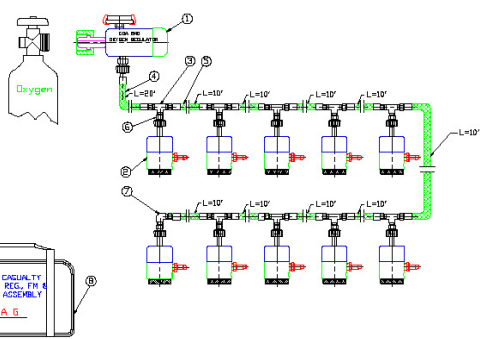

**MCA- 15JM5-20910**  
**MCA- 15JM5-20K10 (with Swivel T's)**  
 10 Flow Meters (0-15 LPM) @ 10ft Intervals  
 CGA-540 Pressure Reducing Regulator  
 20ft Supply Feed Hose and Carry / Storage Bag

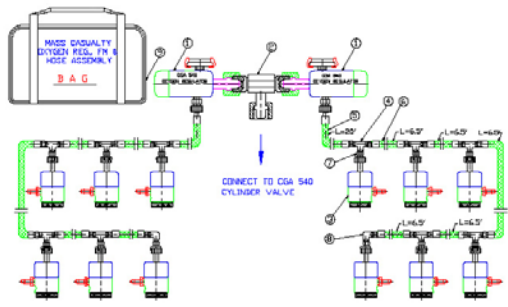

**MCA-75LP7-2AEZ0**  
**MCA-75LP7-2ALZ0 (with Swivel T's)**  
 12 Flow Meters (0-25 LPM) @ 6.5ft Intervals  
 2 CGA-540 Pressure Reducing Regulators  
 20ft Supply Feed Hose and Carry / Storage Bag  
 CGA-540 "T" Adapter

Operational Temperature Range (-25 F to to +150 F)

All MCA Series Multiple Oxygen Flowmeter Distribution Systems can be customized to meet your specific departmental requirements

